



# MANUAL

Standard, lightweight, multifunctional wheelchairs

# MANUAL WHEELCHAIRS



**SPECIFICATIONS**  
AT A GLANCE STARTING ON PAGE 22

## STANDARD FOLDING WHEELCHAIRS

**MEYRA**



STARTING ON PAGE 8  
**SERVICE  
STOCK VERSION**  
3.600



STARTING ON PAGE 9  
**BUDGET  
STOCK VERSION**  
9.050

## LIGHTWEIGHT WHEELCHAIRS

**MEYRA**



PAGE 10  
**HM303D  
STOCK VERSION**  
7.010



PAGE 11  
**HM03MS  
STOCK VERSION**  
7.010



PAGE 12  
**FORMAT  
STOCK VERSION**  
3.940



PAGE 13  
**EUROCHAIR²  
STOCK VERSION**  
2.750



PAGE 16  
**EUROCHAIR² PRO  
STOCK VERSION**  
2.850



PAGE 17  
**EUROCHAIR  
STOCK VERSION**  
1.750

## MULTIFUNCTIONAL WHEELCHAIRS

**MEYRA**



PAGE 18  
**SOLERO LIGHT  
STOCK VERSION**  
9.072

## XXL WHEELCHAIRS

**MEYRA**



PAGE 20  
**EUROCHAIR² XXL  
STOCK VERSION**  
2.850



PAGE 21  
**EUROCHAIR² HD  
STOCK VERSION**  
2.850

THE MEYRA SEAL OF  
INDIVIDUALITY




# MY INDIVIDUAL

**Individual solutions for special requirements** – In our efforts to ensure mobility for each individual, customised product solutions have always been a top priority at MEYRA. As individual as the requirements are, so too are the solutions we develop. Whether for particularly heavy users with a body weight of up to 300 kg or users with specific limitations of their musculoskeletal system – we create the solution that fits, tailored entirely to the handicap in question.

With a highly qualified team that is exclusively responsible for individual solutions, we have expanded the INDIVIDUAL range even further over the course of the last few years. Several thousand individual solutions have thus been created alongside series production. Each request is planned in detail in cooperation with therapists and orthopaedic technicians, and professionally put into effect by our manufacturing team consisting of engineers and technicians. Manufactured with modern CAD equipment – for the accustomed MEYRA quality.

## ORDERING MADE SIMPLE:

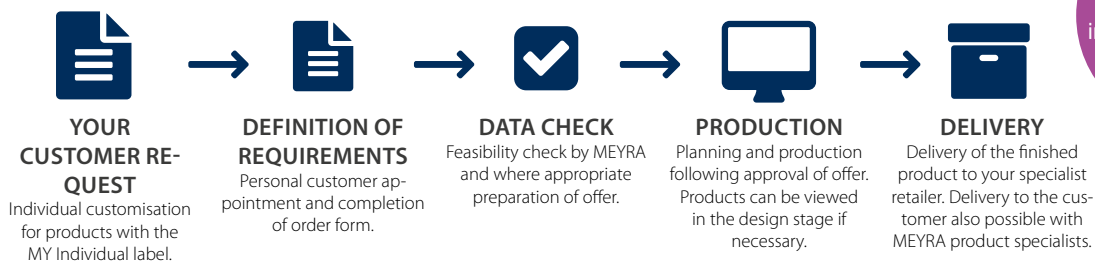

SUBMIT A  
REQUEST:  
[info@meyragroup.com](mailto:info@meyragroup.com)

# GUARANTEE PERIODS & ISO CERTIFICATION

## GUARANTEE TERMS FOR MANUAL WHEELCHAIRS

Precondition for the guarantee on frame and scissor mechanism is the correct use of the wheelchair, which implies regular, professional maintenance. In particular, damage caused by loosened screws or nuts as well as worn fastening holes for the seat belt and back strap due to frequent changing of covers shall be precluded. Furthermore, no changes may have been made to load-bearing parts.

The guarantee excludes damage to the surface, in particular scratches in the paint or other mechanical surface damage and contamination.



## NOTE

It is important to distinguish clearly between the terms "permissible user weight" and "permissible total weight" of a wheelchair, the latter being the sum of:

- a) permissible user weight
- b) tare weight of the wheelchair
- c) additional load

MODEL TYPE	MODEL NO.	PERMISSIBLE USER WEIGHT	ISO 10542-2	ISO 7176-19	GUARANTEE ON FRAME & SCISSOR MECHANISM*
<b>SERVICE Stock version</b>	3.600	120 kg	-	-	3 years
<b>BUDGET Stock version</b>	9.050	130 kg	-	✓	3 years
<b>HM303D Stock version</b>	7.010	110 kg	-	-	1 year
<b>HM03MS Stock version</b>	7.010	100 kg	-	-	1 year
<b>FORMAT Stock version</b>	3.940	130 kg	-	✓	3 years
<b>EUROCHAIR<sup>2</sup> Stock version</b>	2.750	130 kg	-	✓	3 years
<b>EUROCHAIR<sup>2</sup> PRO</b>	2.850	130 kg	-	✓	4 years
<b>EUROCHAIR<sup>2</sup> XXL</b>	2.850	160 kg / 200 kg	-	✓	3 years
<b>EUROCHAIR<sup>2</sup> HD</b>	2.850	300 kg	-	-	-
<b>EUROCHAIR stock version</b>	1.750	130 kg	-	✓	3 years
<b>SOLERO LIGHT Stock version</b>	9.072	130 kg	-	✓	3 years

## MEYRA QUALITY

### TEST TRACKS & CRASH TESTS

All wheelchairs are fully tested on a test track. A wheelchair is only approved for delivery having successfully completed the internal quality tests. In addition, a crash test is carried out on new models before they are launched on the market to ensure that the wheelchair can be used in place of a car seat. Go to our website to find out more about crash tests: [www.meyra.com/information-center/laws-and-regulations/crash-test](http://www.meyra.com/information-center/laws-and-regulations/crash-test)

### THE MEYRA ENDURANCE TEST

MEYRA's tests go beyond the level of standards. We carry out stress and reliability tests with increased loads over and above the level required by standards. This is unique and ensures that our wheelchairs meet our high quality requirements.

REQUIRED BY LAW:	OUR TESTING STANDARD:
Double-drum test according to legal standard DIN EN 21856 DIN EN 12183	Double-drum test according to legal standard and additional MEYRA test procedure ISO 7176-8
<b>200,000 cycles</b>	<b>200,000 – 300,000 cycles</b>
<b>without overload</b>	<b>with 10% overload</b>
-	<b>additional stress test</b>





# MEDICAL CONDITIONS & FUNCTIONAL REQUIREMENTS



## MEDICAL CONDITIONS

Traumatic brain injury MS, ALS, severe dementia	-	-	-	-	-	-	-	-	-	-	-	-
Stroke Parkinson's	-	-	-	-	-	-	-	-	-	-	-	-
Leg amputations with degenerative diseases	-	-	-	-	-	-	-	-	-	-	-	-
Geriatric disorders long-term fittings	-	-	-	-	-	-	-	-	-	-	-	-
Obesity with sequelae	-	-	-	-	-	-	-	-	-	-	-	-
Post-OP conditions or traumas short-term fittings	-	-	-	-	-	-	-	-	-	-	-	-
Geriatric disorders short-term fittings	-	-	-	-	-	-	-	-	-	-	-	-
Clinic and institutional use for transfer and short-term fittings	-	-	-	-	-	-	-	-	-	-	-	-
	SERVICE 3.600 Stock version	BUDGET 9.050 Stock version	FORMAT 3.940 Stock version	EUROCHAIR 1.750 Stock version	HM303D 7.010 Stock version	HM03MS 7.010 Stock version	EUROCHAIR² 2.750 Stock version	EUROCHAIR² PRO 2.850	EUROCHAIR² XXL 2.850	EUROCHAIR² HD 2.850	SOLERO LIGHT 9.072 Stock version	
FUNCTIONAL REQUIRE- MENTS IN THE VARIOUS MEDICAL FIELDS	Standard	Standard lightweight	Lightweight Stock version				Lightweight , individual + XXL				Multi- function	
Palliative care	-	-	-	-	-	-	-	-	-	-	-	-
Geriatrics	-	-	-	-	-	-	-	-	-	-	-	-
Neurology	-	-	-	-	-	-	-	-	-	-	-	-
Orthopaedics, traumatology	-	-	-	-	-	-	-	-	-	-	-	-
Internal medicine	-	-	-	-	-	-	-	-	-	-	-	-



## FIND OUT MORE

AT [WWW.MEYRA.COM](http://WWW.MEYRA.COM) OR CONTACT  
YOUR LOCAL SPECIALIST RETAILER


# COMMON PARTS







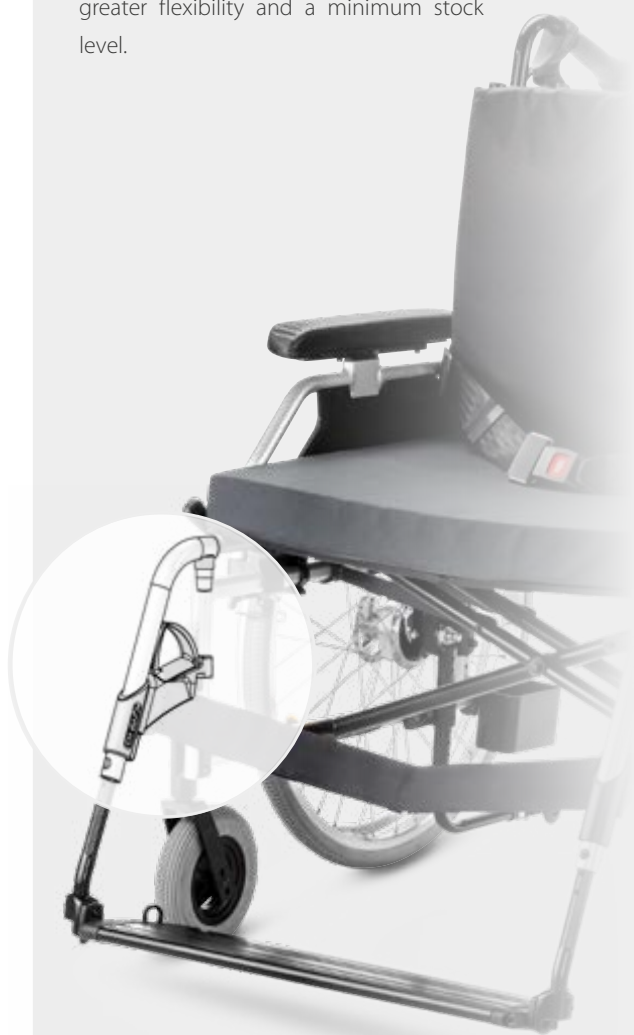

## PARTS COMPATIBILITY

Very economical due to maximum compatibility of parts – the idea behind MEY-RA's modular system. The modular system is already an integral feature in the EURO-CHAIR family and has also established itself in the rest of the manual sector. The widespread uniformity of the individual parts used in our wheelchairs also guarantees greater flexibility and a minimum stock level.

## COMPONENTS

Side panel, standard
Side panel, height-adj.
Pressure brake
Brake lever extension
Drum brake <sup>2)</sup>
Footrest
Swivel-up footrest
Amputee leg support
Back tube
Back upholstery
Seat upholstery <sup>1)</sup>
Cross brace
Castor wheels
Drive wheels
Support castor wheels
Therapy table
Therapy table, half lap
Walking aid holder
Retainer strap
Passive lighting
Spoke guards



							
SERVICE Stock version 3.600	BUDGET Stock version 9.050	FORMAT Stock version 3.940	EUROCHAIR Stock version 1.750	EUROCHAIR² Stock version 2.750	EUROCHAIR² PRO 2.850	EUROCHAIR² XXL 2.850	SOLERO LIGHT Stock version 9.072
-	○	○	○	○	○	○	-
-	○	○	○	○	○	○	-
○	○	○	○	○	○	○	-
○	○	○	○	○	○	○	-
-	○	○	○	○	○	○	-
-	○	○	○	○	○	○	○
-	○	○	○	○	○	○	○
-	○	○	○	○	○	○	○
-	○	○	○	○	○	○	○
-	-	-	-	○	○	○	-
-	-	-	○	○	○	○	-
-	-	-	○	○	○	○	-
-	-	-	○	○	○	○	-
○	○	○	○	○	○	○	○
○	○	○	○	○	○	○	○
-	○	○	○	○	○	○	-
○	○	○	○	○	○	up to SW 53 cm	○
○	○	○	○	○	○	○	○
-	○	○	○	○	○	○	-
○	○	○	○	○	○	○	○
-	○	○	○	○	○	○	-
○	○	○	○	○	○	○	○

<sup>1)</sup> Models 1.XXX use thinner fastening screw

<sup>2)</sup> poss. different mount

# SERVICE STOCK VERSION 3.600

## USER PROFILE

- ✓ Geriatric fitting
- ✓ Rheumatism
- ✓ As transport or pushing aid in hospitals, clinics and care homes
- ✓ For short-term use around the home or following an accident
- ✓ Transport of guests with restricted mobility, e.g. in hotels
- ✓ As a walking aid (walking while seated in wheelchair)

## FOR UNIVERSAL USE

The SERVICE is ideally suited for non-permanent use (e.g. following an operation) as a transport and mobility aid.

- Possibilities for adaptation to the body measurements of the user
- Extensive range of accessories
- Smooth folding due to scissor guide tubes, thus high-quality folding system with long service life and good parallel guidance of the side frames
- Hand rim propulsion/push mode
- Side panels, armrests and legrests detachable
- Knee lever brake
- Clinic version available in chrome
- PU armrests suitable for therapy table



Castor wheel puncture-proof 200 x 50 mm



Detachable, swing-away legrests



Side panel chamfered at the front with padded armrest




Aluminium hand rim with 15 mm clearance

Seat height front / rear: 520 mm / 490 mm  
Seat depth: 440 mm

Frame colour Chrome or Silverline



User weight up to 120 kg

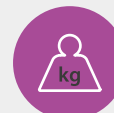

Empty weight of wheelchair from 18.5 kg up



Seat widths from 400 – 510 mm





# BUDGET STOCK VERSION 9.050

## USER PROFILE

- ✓ Push-hand rim drive
- ✓ Rheumatism
- ✓ Geriatric fitting.

## ECONOMICAL AND INNOVATIVE

Thanks to its ease of use in everyday situations, the BUDGET is well suited for use around the home and in inpatient care.

- Front seat height variable without swapping parts
- Comfort-remove side panel
- Standard wheelbase extension with unique optimisation of manoeuvrability for attendant
- Folding
- Hand rim propulsion/push mode
- Legrests swing-away inwards and outwards, angle-adjustable footplates
- Lockable, backward-folding armrests
- All drive wheels incl. quick release axle, stepper aid as standard
- Pressure or drum brake for attendant (option)



Knee lever brake



Lockable,  
backward-folding  
armrests



Detachable, swing-away  
legrests



Back with lumbar contour



Comfort-remove side panel



Front seat height-adjustable without replacing parts  
Seat height front / rear:  
470 – 510 mm / 440 – 490 mm  
Seat depth: 430 mm



Frame colour Silverline



User  
weight up to  
130 kg



Empty weight  
of wheelchair  
from 19.6 kg up



Seat widths  
from  
380 – 510 mm



# HM303D STOCK VERSION 7.010

## USER PROFILE

- ✓ As transport or pushing aid in hospitals, clinics and care homes
- ✓ For short-term use around the home
- ✓ Transport of guests with restricted mobility, e.g. in hotels

## COMPACT AND EASY TO CARRY – PRODUCT WEIGHT FULLY ASSEMBLED: ONLY 10 KG

The product is a lightweight, folding wheelchair weighing only 10 kg that offers additional safety thanks to its dual brake control. It has 22-inch wheels with quick-release axles, which can be detached to facilitate transfer and make storage and transport more convenient. The swing-away and detachable legrests offer more freedom of movement. The lightweight design of the wheelchair allows the user to carry it easily and lift it into a car with minimal effort. Its folding system and hinged backrest make it easy to transport. It is ideal for out and about and for users seeking independence and flexibility.

- Frame colour: matt black
- Portable
- Dual brake system
- Transport wheelchair



Detachable, swing-away legrests



Dual brake system



Back bag

Folding back

22" drive wheels



User weight up to 110 kg



Empty weight of wheelchair only 10 kg



Seat width 450 mm

# HM03MS STOCK VERSION 7.010


## USER PROFILE

- ✓ As transport or pushing aid in hospitals, clinics and care homes
- ✓ For short-term use around the home
- ✓ Transport of guests with restricted mobility, e.g. in hotels


## THE SUPER-LIGHTWEIGHT, FOLDING TRANSPORT WHEELCHAIR WEIGHS ONLY 9.6 KG

Weighing in at only 9.6 kg, this wheelchair is lightweight and can be folded, making it ideal for transporting. It has small 12" drive wheels and also offers a dual brake system for the user and attendant. It is specially designed to meet the individual needs of specific users. Swivel-up legrests and armrests and a folding backrest make it compact to pack. The wheelchair offers easy manoeuvrability and stability. Perfect for anyone looking for mobility and flexibility in a compact format.

- Frame colour: navy blue
- Portable
- Dual brake system
- Transport wheelchair

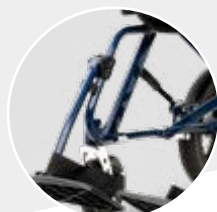

Swivel-up  
legrests

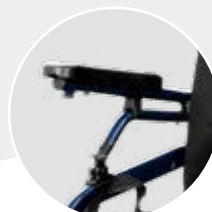

Swivel-up  
armrests



Strong cross brace



Folding back

12" drive wheels



User  
weight up to  
100 kg



Empty weight  
of wheelchair  
only 9.6 kg



Seat width  
450 mm

# FORMAT STOCK VERSION 3.940

## USER PROFILE

- ✓ At home and in inpatient care
- ✓ Geriatric fitting
- ✓ Rheumatism

## THE QUALITY ALTERNATIVE FOR REUSE

The FORMAT is perfect for getting started with independent mobility. If the user requires additional adjustment, the FORMAT provides support when standing up and transferring, as the seat height can be adjusted to the user's height.

- Folding
- High-quality adjustable back
- Height-adjustable push handles
- Height-adjustable side panels
- Legrests swing-away inwards and outwards
- Adaptable seat height
- Height-adjustable, lockable, backward-folding armrests, tuck in behind the back
- All drive wheels incl. quick release axle
- Knee lever brake or drum brake for attendant (option)
- Angle-adjustable footplate, stepper aid and wheelbase extension as standard
- Links in with MEYRA's EUROCHAIR Family modular system



Height-adjustable  
push handles as standard



Height-adjustable  
side panels as standard



Safety package (optional): passive lighting, support castors, safety belt, Walking stick holder



Height-adjustable push handles as standard

High-quality adjustable back as standard

Height-adjustable side panels as standard

Seat height front / rear:  
470 – 510 mm / 440 – 490 mm  
Seat depth: 400 – 430 mm

24" aluminium drive wheel

Detachable, swing-away legrests



User  
weight up to  
130 kg



Empty weight  
of wheelchair  
from 17.5 kg up



Seat widths  
from  
380 – 510 mm





# EUROCHAIR<sup>2</sup> STOCK VERSION 2.750

## THE LIGHTWEIGHT CHAIR FOR UNIVERSAL USE

The EUROCHAIR<sup>2</sup> encourages independent mobility by means of fine adjustment of the seat, backrest and chassis to the needs of the user. Ease of use, improved smoothness with low propelling force, and reduced weight facilitate everyday handling.

- Folding and individually adaptable
- Side panel height-adjustable, swivel-up, swing-away, detachable
- Angle-adjustable steering head
- Height-adjustable side panel with depth-adjustable armrest
- Adjustable back
- MEYRA modular system ensures extensive range of accessories and quick supply of spare parts

### USER PROFILE


- ✓ For permanent everyday use
- ✓ For optimum fitting for a wide range of medical conditions through high adjustment variability
- ✓ Primarily for use at home, but also in rehabilitation clinics for achieving therapeutic goals
- ✓ For taller and larger users



Safety package and  
seat cushion included



Height-adjustable  
side panel with  
depth-adjustable  
armrest

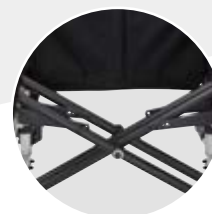

Double scissor as  
standard



Seat belt and back strap: Padded, black

Seat height variable: 470 – 520 mm  
Seat depth variable: 400 – 490 mm

Wheelbase extension by turning the  
varioblock backwards

Detachable, swing-away legrests. Footplates  
angle-adjustable, incl. single calf strap



User  
weight up to  
130 kg



Empty weight  
of wheelchair  
from 17.5 kg up



Seat widths  
from  
380 – 530 mm











HIGH-QUALITY.  
DURABLE.  
MODULAR.

 MEYRA®

# EUROCHAIR<sup>2</sup> PRO 2.850

## USER PROFILE

- ✓ Fully equipped for permanent use
- ✓ Optimum fitting for a wide range of medical conditions owing to high adjustment variability
- ✓ Primarily for use at home, but also in rehabilitation clinics for achieving therapeutic goals
- ✓ For use with various drive types
- ✓ HEMI components

## FUNCTIONAL AND INTELLIGENT

The EUROCHAIR<sup>2</sup> PRO encourages independent mobility by means of fine adjustment to the user's needs. It offers low-fatigue sitting over long periods, easy operation of side panels and legrests, and smooth running properties requiring minimal propelling force.

- Side panel height-adjustable, swivel-up, swing-away, detachable
- Double scissor as standard
- Folding
- Legrests swing-away inwards and outwards
- Armrests detachable, backward folding
- Knee lever brake or drum brake for attendant (option)
- Stepper aid as standard
- Easy wheelbase extension by turning the varioblock
- Very well suited to various areas of use, as the same seat height can be set with different wheel sizes



Castor wheel 180 mm PU



Height-adjustable side panel with depth-adjustable armrest



Wheelbase extension by turning the varioblock backwards



Adjustable back standard



Seat depths from 400 – 490 mm  
Seat height variable



Detachable, swing-away legrests. Footplates angle-adjustable, incl. single calf strap



User weight up to 130 kg



Empty weight of wheelchair from 17.5 kg up



Seat widths from 300 – 530 mm



# EUROCHAIR STOCK VERSION 1.750

## USER PROFILE

- ✓ For permanent everyday use
- ✓ For optimum fitting for a wide range of medical conditions through high adjustment variability
- ✓ Primarily for use at home, but also in rehabilitation clinics for achieving therapeutic goals
- ✓ For taller and larger users

## INDIVIDUALLY ADAPTED

The EUROCHAIR presents numerous possibilities for adjustment without having to swap parts and is therefore ideal for individual adaptation. The position of the drive wheel can be adjusted horizontally and vertically for even greater flexibility. Very short modification times, high variability thanks to retrofit kits, and low weight make everyday handling easier.

- Stock version
- Optimal fitting for numerous medical conditions
- High degree of adjustment variance and individual adaptation
- Service-friendly
- Numerous possibilities for adjustment without swapping parts
- Easy to transport
- Extensive modular system
- Folding



Angle-adjustable  
steering head, aluminium  
Castor fork



High variability thanks to  
retrofit kits



Fast and precise conver-  
sion without swapping  
parts



Flip-up, swing-away and detachable side panel

Seat height front / rear:  
430 – 510 mm / 430 – 530 mm  
Seat depth: 400 – 460 mm

Two frame colours to choose from

Legrests swing-away inwards and outwards,  
angle-adjustable footplates



User  
weight up to  
130 kg



Empty weight  
of wheelchair  
from 15 kg up  
(with pressure  
brake)



Seat widths  
from  
380 – 530 mm





# SOLERO LIGHT STOCK VERSION 9.072

## USER PROFILE

- ✓ Geriatric fitting in retirement and care centres as well as at home
- ✓ For physical dysfunctions that require a good, supportive sitting position
- ✓ Good suitability and mobilisation for self-propelling wheelchair users

## EASY HANDLING

The very narrow overall width of the SOLERO LIGHT is beneficial when passing through tight doorways or in confined spaces.

- Simple to operate
- Infinitely variable adjustment of the seat and back angle
- Pressure brake, optional drum brake
- Assembly of all auxiliary drives possible
- Gas pressure spring ensures gentle recline and pull-up
- Compact and low seat height
- Variable seat width and seat depth adjustment
- Headrest height, depth and angle-adjustable
- Armrests height-adjustable and detachable with memory function
- Centre of gravity oriented frame design provides increased structural stability



Armrest height and depth-adjustable



Transverse push bar (optional)



Headrest with adjustable side padding and ball joint (optional)



Seat tilt from -2 to +25°

Seat height with seat padding: 440 / 520 mm  
Seat depth: 420 – 480 mm

Legrests angle-adjustable and lockable, detachable, swing-away and swivel-up



User weight  
130 kg



Empty weight  
of wheelchair  
from 38 kg up



Seat widths  
from  
430 – 460 mm /  
480 – 510 mm







# XXL FOR SPECIAL REQUIREMENTS

MORE XXL PRODUCTS  
AT [WWW.MEYRA.COM](http://WWW.MEYRA.COM) OR IN  
OUR BARIATRICS BROCHURE.



# EUROCHAIR<sup>2</sup> XXL 2.850

## USER PROFILE

- ✓ Suitable for heavy users up to 200 kg
- ✓ Long-term fitting with high adaptation requirement for independent mobility
- ✓ Primarily for use at home, but also in rehabilitation clinics for achieving therapeutic goals
- ✓ Use with various drive types

## MAXIMUM STABILITY

The EUROCHAIR<sup>2</sup> XXL makes driving easier thanks to numerous adjustment options on the chassis. An ergonomic gripping angle relieves the strain on the shoulder joints. Improved efficiency increases independence and supports optimum utilisation of available propelling force. This leads to easier mobility and greater manoeuvrability.

The addition of a transverse pushbar and/or a one-piece footboard gives the user the greatest possible support in the sitting position.

- Folding
- Particularly stable and available in five seat widths
- Adjustable back
- Drum brake for attendant (option)
- Stepper aid as standard



One-piece footboard



Increased stability with the use of a double scissor mechanism



Height-adjustable side panel with depth-adjustable armrests as standard



Seat height front / rear variable:  
460 – 520 mm  
Seat depth: 430 / 460 / 490 mm

Extremely sturdy

Detachable, swing-away legrests



User weight up to 160 kg (reinforced version up to 200 kg)



Empty weight of wheelchair from 19 kg up



Seat widths from 480 – 650 mm



# EUROCHAIR<sup>2</sup> HD 2.850

## SPECIAL REQUIREMENTS

The EUROCHAIR<sup>2</sup> HD enables a heavy person to move independently and/or be cared for and transported in barrier-free indoor and outdoor areas. Good adaptation of the seat, backrest and chassis to the user supports independent mobility as far as possible. The transverse pushbar and the one-piece footboard give the user the greatest possible support in the sitting position.


### USER PROFILE

- ✓ Suitable for heavy users up to 300 kg
- ✓ Facilitates care and transfer for the user and the attendant
- ✓ Use in retirement and nursing homes, therapy and rehabilitation centres or around the home to support therapeutic goals.


- Reinforced back, seat belt and rear side panel mount
- Folding
- High-strength steel frame
- One-piece footboard
- Armrests detachable, backward folding
- Drum brake for attendant
- Stepper aid as standard



One-piece footboard

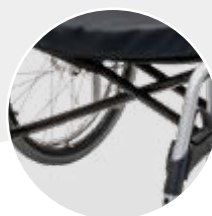

Increased stability  
with the use of a  
double scissor mechanism

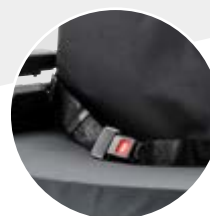

Reinforced seat belt



Seat height front / rear:  
500 / 475 mm or 475 / 450 mm  
Seat depth: 460 / 500 mm

Reinforced version up to 300 kg user weight  
included in scope of delivery:

- Reinforced back strap
- Reinforced seat belt
- One-piece footboard
- Drum brake
- Transverse pushbar



User  
weight up to  
300 kg



Empty weight  
of wheelchair  
from 35 kg up



Seat widths  
from  
700 – 800 mm




# SPECIFICATIONS

Dimensional tolerances  $\pm 10$  mm,  $\pm 2^\circ$ , data subject to design changes

	SERVICE STOCK VERSION 3.600	BUDGET STOCK VERSION 9.050	HM303D STOCK VERSION 7.010	HM03MS STOCK VERSION 7.010	FORMAT STOCK VERSION 3.940
Seat width in mm	400 / 450 <sup>3)</sup> / 480 / 510	380 / 400 / 430 / 460 / 480 / 510	450	450	380 / 400 / 430 / 460 / 480 / 510
Width, ready to move, in mm	620 / 670 / 700 / 730	560 / 580 / 610 / 640 / 660 / 690	645	600	560 / 580 / 610 / 640 / 660 / 690
Width, folded, in mm	310	280	-	-	280
Seat depth in mm	440	430	410	410	400 – 430
Seat height in mm	-	-	450	465	-
Seat height, front, in mm	520	470 – 510	-	-	470 – 510
Seat height, rear, in mm	490	440 – 490	-	-	440 – 490
Total height in mm	920	960	860	875	960
Back height in mm	400	420	400	410	420
Back angle in degrees	-	-	$82^\circ \pm 1^\circ$	$82^\circ \pm 1^\circ$	-
Seat tilt in degrees	-	-	-	-	-
Armrest height in mm	240	210.5	220	220	215
Footboard length in mm	-	-	390	390	-
Length with / without footrest	1,080 / 780	1,040	942 / -	1,075 / -	1,040
User weight in kg	120	130	110	100	130
Empty weight in kg	19	18.5	10	9.5	17.5
Permissible total weight in kg	139	148.5	-	-	147.5
Transport weight in kg	11	11.5	-	-	10.5
Drive wheel in mm (inches)	610 (24")	610 (24")	22"	12"	610 (24")
Static stability (also with brake applied) in degrees	-	-	-	-	-
Brake test according to STANDARD in degrees	-	-	-	-	-
Pack size in mm	-	-	930 x 360 x 710	910 x 350 x 720	-

EUROCHAIR <sup>2</sup> STOCK VERSION 2.750	EUROCHAIR <sup>2</sup> PRO 2.850	EUROCHAIR STOCK VERSION 1.750	SOLERO LIGHT STOCK VERSION 9.072	EUROCHAIR <sup>2</sup> XXL 2.850	EUROCHAIR <sup>2</sup> HD 2.850
380 / 400 / 430 / 460 / 480 / 500 / 530	300 / 320 / 340 / 360 / 380 / 400 / 430 / 460 / 480 / 500 / 530	380 / 400 / 430 / 460 / 480 / 500 / 530	430 – 460 / 480 – 510	480 / 500 / 530 / 580 / 650	700 / 750 / 800
580 / 600 / 630 / 660 / 680 / 700 / 730	500 / 520 / 540 / 560 / 580 / 600 / 630 / 660 / 680 / 700 / 730	560 / 580 / 610 / 640 / 660 / 680 / 710	min. 620	680 / 700 / 730 / 780 / 850	920 / 970 / 1,020
310	310	280		310	310
variable 400 – 490	400 / 430 / 460 / 490	400 – 460	420 – 480	430 / 460 / 490	460 / 500
-	-	-	440 / 520	-	-
variable 470 – 520	variable 390 – 560	variable 430 – 510	-	variable 460 – 520	475 / 500
	variable 390 – 520	variable 430 – 530	-	variable 460 – 520	475 / 450
930	variable 790 – 1,050	960	1,040 – 1,180	930	960
400 / 420 / 440	400 / 420 / 440 / 460 / 480 / 500	400 / 420	530 – 560	400 / 420 / 440 / 460 / 480 / 500	430
-	-	-	0° – 45°	-	-
-	-	-	(-2°) – 25°	-	-
210 – 260	210 – 260	from 230	210 – 270	210 – 260	230 – 330
-	-	-	-	-	-
1,050 <sup>4)</sup>	1,050	from 1,000	1,180 / 910	max. 1,070 / max. 840	max. 1,140 / max. 880
130	130	130	130	160 / 200 / 300	300
17.5	17.5	15 <sup>5)</sup>	38	19	34
147.5	147.5	150	168	179 / 219	334
-	10	-	23	12	27.5
616 (24")	22" 24"	600 (24")	610 (24")	610 (24")	610 (24")
-	-	-	-	-	-
7°	-	7°	-	-	-
-	-	-	-	-	-





*inspire  
joy of life*

Issue 2025/01

**MEYRA GmbH** Kleiststraße 49 | 32457 Porta Westfalica | Germany  
info@meyragroup.com | Tel.: +49 5733 922 - 0 | Fax: +49 5733 922 - 9311 | [www.meyra.com](http://www.meyra.com)

2025-08 We reserve the right to make technical modifications to our products  
and assume no liability for printing errors or variation in colour in our printed matter.

Ident-Nr. 261 505 601

Picture credits: Page 1,14,15,24 © "Halfpoint" – stock.adobe.com; Page 3 © "Gorodenkoff" – stock.adobe.com;  
Page 19 © "FredFroese" – istockphoto.com